Simba

Gebrauchsanweisung
Direction for Use
Notice d'utilisation
Manuale d'uso
Gebruikershandleiding
Manuale Per l'Uso
Manual do Utilizador
Brukerveiledning
Bruksanvisning
Brugervejledning

zippie®

Vorwort

Liebe Kundin/Lieber Kunde,

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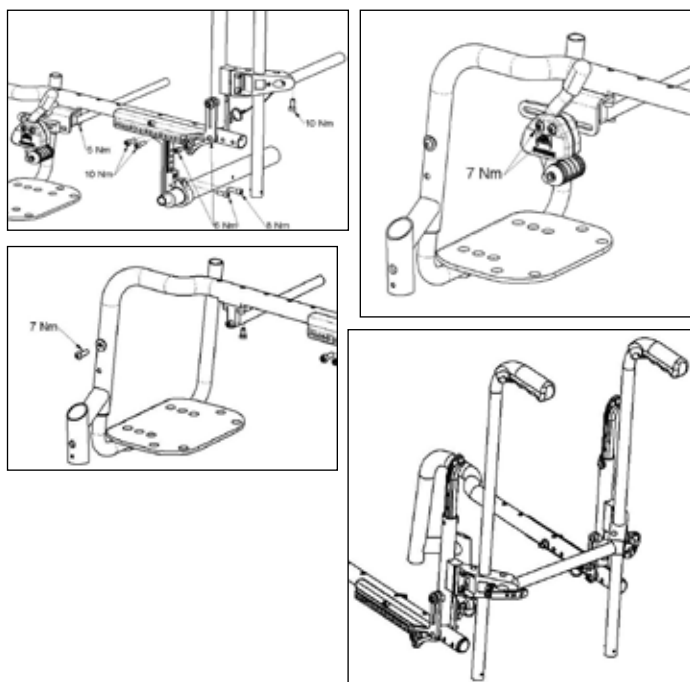
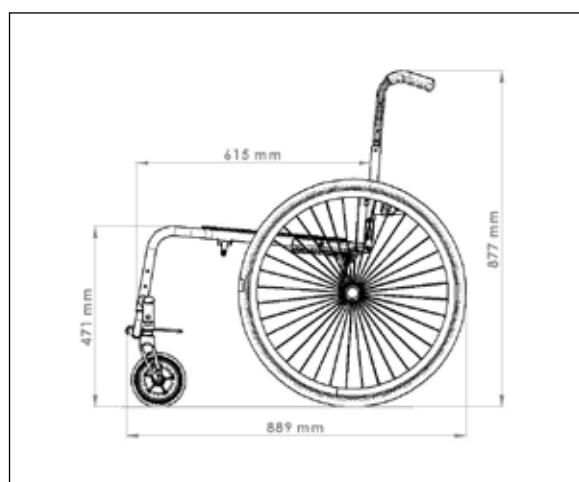
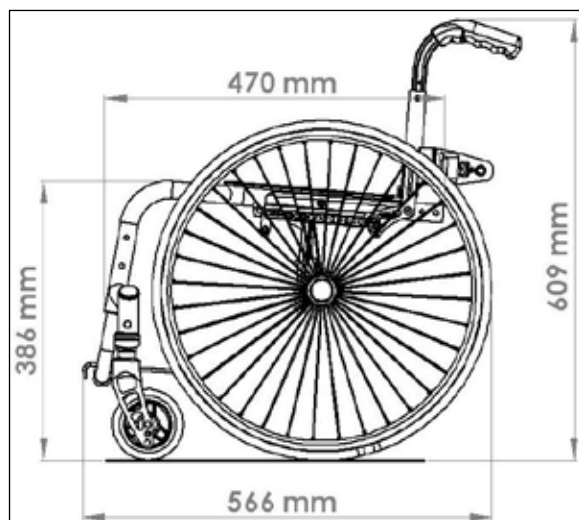
Foreword

Dear Customer,

We are very pleased that you have decided in favour of a high-quality SUNRISE MEDICAL product.

This user's manual will provide numerous tips and ideas so that your new wheelchair can become a trustworthy and reliable partner in your life.

For Sunrise Medical, it is very important that we have a good relationship with our customers. We like to keep you up-to-date about new and current developments at our company. Keeping close to our customers means: fast service, as little red tape as possible, working closely with customers. When you need replacement parts or accessories, or if you just have a question about your wheelchair – we are there for you.

We want you to be satisfied with our products and service. At Sunrise Medical we are constantly working to develop our products further. For this reason, changes can occur in our palette of products with regard to form, technology and equipment. Consequently, no claims can be construed from the data or pictures contained in this user's manual.

The management system of SUNRISE MEDICAL is certified to DIN ISO 9001, ISO 13485 and ISO 14001.

As the manufacturer, SUNRISE MEDICAL, declares that the lightweight wheelchairs conform to the 93/42/EEC / 2007/47/EEC guideline.

Please contact your local, authorised SUNRISE MEDICAL dealer if you have any questions regarding the use, maintenance or safety of your wheelchair.

In case there is no authorised dealer in your area or you have any questions, contact Sunrise Medical either in writing or by telephone.

Sunrise Medical LTD.
Sunrise Business Park
High Street, Wollaston
West Midlands DY8 4PS
England

Phone: +44 (0) 1384 44 66 88
Fax: +44 (0) 1384 44 66 99

Sunrise Medical is ISO 9001 certified, which ensures quality at all stages of the development and production of this wheelchair.

This product is manufactured to comply with Medical Device Directive 93/42/EEC.

Your wheelchair is not suitable for use as a seat in vehicles. Do not sit in the wheelchair during transportation. Use the vehicle seat or a properly secured transport approved seat.

IMPORTANT:
DO NOT USE YOUR WHEELCHAIR UNTIL THIS MANUAL HAS BEEN READ AND UNDERSTOOD.

Table of contents

Table of contents	15
Foreword about wheelchairs	15
General Safety Notes and Driving Limits	16 - 17
Guarantee	17
Wheelchair Components	18
Handling	19
Step tubes	19
Brakes	19
Footplate adjustment	19
Castors	20
Seat Depth	20
Wheel Alignment	20
Backrest	21
Armrest	21
Sideguard	21
Push Handles	21
Positioning belt	21
Anti-tip Tubes	22
Seat	22
Crutch Holder	22
Tyres and Mounting	22
Maintenance and Care	22
Lap belt	23
Trouble-shooting	24
Disposal / Recycling of Materials	24
Nameplate	24
Torque	25
Technical Data	25 - 26

Use

The children's wheelchair is exclusively for a child who is unable to walk or has limited mobility, for their own personal use in- and outdoor.

The maximum weight limit (includes the user and any weight of accessories fitted to the wheelchair) is marked on the serial number label, which is affixed to the crossbar or stabiliser bar below the seat.

Warranty can only be taken on if the product is used under the specified conditions and for the intended purposes.

The expected life of the wheelchair is 5 years.

Please DO NOT use or fit any 3rd party components to the wheelchair unless they are officially approved by Sunrise Medical.

Area of application

The variety of fitting variants as well as the modular design mean that it can be used by those who cannot walk or have limited mobility because of:

- Paralysis
- Loss of extremity (leg amputation)
- Extremity defect deformity
- Joint contractures/joint injuries
- Illnesses such as heart and circulation deficiencies, disturbance of equilibrium or cachexia.

When considering provision, please also note the body size, weight, physical and psychological constitution, the age of the person, living conditions and environment.

General safety notes and driving limits

The engineering and construction of this wheelchair has been designed to provide maximum safety. International safety standards currently in force have either been fulfilled or exceeded in parts. Nevertheless, users may put themselves at risk by improperly using their wheelchairs. For your own safety, the following rules must be strictly observed.

Unprofessional or erroneous changes or adjustments increase the risk of accident. As a wheelchair user, you are also part of the daily traffic on streets and pavements, just like anyone else. We would like to remind you that you are therefore also subject to any and all traffic laws.

Be careful during your first ride in this wheelchair. Get to know your wheelchair.

Before each use, the following should be checked:

- Quick-release axles on the rear wheels
- Velcro on seats and backrests
- Tyres, tyre pressure and wheel locks.

Before changing any of the adjustments of this wheelchair, it is important to read the corresponding section of the user's manual.

It is possible that potholes or uneven ground could cause this wheelchair to tip over, especially when riding uphill or downhill.

When riding over a step or up an incline frontally, the body should be leaning forward.

Danger!

NEVER exceed the maximum load of 125 kg for driver plus any items carried on the wheelchair. If you exceed the maximum load, this can lead to damage to the chair, or you may fall or tip over, lose control and may lead to serious injury of the user and other people.

Danger!


When it is dark, please wear light clothing or clothing with reflectors, so that you can be seen more easily. Make sure that the reflectors on the side and back of the wheelchair are clearly visible. We would also recommend that you fit an active light.

Danger!

To avoid falls and dangerous situations, you should first practice using your new wheelchair on level ground with good visibility.

Danger!

When getting on or off the wheelchair, do not use the footboards. These should be flipped up beforehand and swung to the outside as far as possible. Always position yourself as close as possible to the place where you wish to sit.

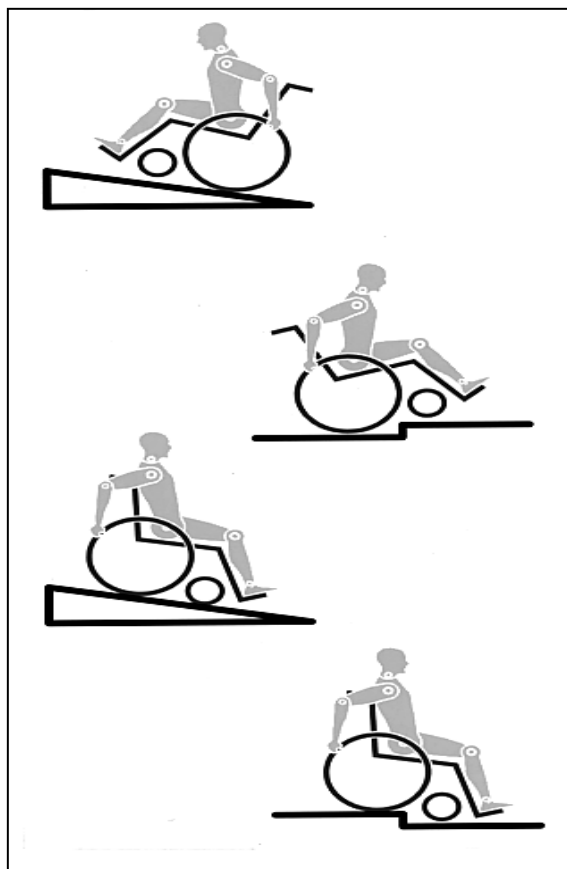Danger!

Only use your wheelchair properly. For example, avoid travelling against an obstacle without braking (step, kerb edge) or dropping down gaps.

Danger!

The wheel locks are not intended to brake your wheelchair. They are only there to ensure that your wheelchair does not begin rolling unintentionally. When you stop on uneven ground, you should always use the wheel locks to prevent such rolling. Always apply both wheel locks; otherwise, your wheelchair could tip over.

Danger!

Explore the effects of changing the centre of gravity on the behaviour of the wheelchair, for example on inclines, slopes, all gradients or when overcoming obstacles. Do this with the secure aid of a helper.

Danger!

With extreme settings (e.g. rear wheels in the most forward position) and less than perfect posture, the wheelchair may tip over even on a level surface.

Danger!

Lean your upper body further forward when going up slopes and steps.

Danger!

Lean your upper body further back when going down slopes and steps. Never try to climb and descend a slope diagonally.

Danger!

Avoid using an escalator which may lead to serious injury in the event of a fall.

Danger!

Do not use the wheelchair on slopes > 10°. The Dynamic safe slope is dependant on the chair configuration, the users abilities and the style of riding. As the users abilities and style of riding cannot be pre-determined then the max safe slope cannot be determined.

Therefore this must be determined by the user with the assistance of an attendant to prevent tipping. It is strongly recommended that inexperienced users have Ant-tips fitted.

Danger!

It is possible that potholes or uneven ground could cause this wheelchair to tip over, especially when riding uphill or downhill.

Danger!

Do not use your wheelchair on muddy or icy ground. Do not use your wheelchair where pedestrians are not allowed.

Danger!

To avoid hand injuries do not grab in between the spokes or between the rear wheel and wheel lock when driving the wheelchair.

Danger!

In particular when using lightweight metal handrims, fingers will easily become hot when braking from a high speed or on long inclines.

Danger!

Only attempt stairs with the help of an attendant. There is equipment available to help you, e.g. climbing ramps or lifts, please use them. If there is no such equipment available, then the wheelchair must be tipped and pushed, never carried, over the steps (2 helpers). We recommend that users over 100 kg in weight do not use this stairway manoeuvre!

In general, any anti-tip tubes fitted must be set beforehand, so that they cannot touch the steps, as otherwise this could lead to a serious tumble. Afterwards the anti-tip tubes must be set back to their correct position.

Danger!

Make sure that the attendant only holds the wheelchair using securely mounted parts (e.g. not on the footrests or the sideguards).

Danger!

When using the lifting ramp make sure that the anti-tip tubes fitted are positioned outside the danger area.

Danger!

Depending on the diameter and setting of the castors, as well as the centre of gravity setting of the wheelchair, the castors may begin to wobble at high speeds. This can lead to the castors being blocked and the wheelchair may tip over. Therefore, please make sure that the castors are adjusted correctly (see the Chapter "Castors"). In particular, do not travel on an incline without brakes, travel at a reduced speed. We recommend that novice users use anti-tip tubes.

Danger!

Anti-tip tubes should prevent the chair tipping over backwards unintentionally. Under no circumstances should they take the place of transit wheels, and be used to transport a person in a wheelchair with the rear wheels removed.

Danger!

When reaching for objects (which are in front of, to the side or behind the wheelchair) make sure that you do not lean too far out of the wheelchair, as if you change the centre of gravity there is a risk of tipping or rolling over. The hanging of additional load (back pack or similar items) onto your chair backposts can affect the rearward stability of your chair, especially when used in combination with recliner backrests. This can cause the chair to tip backwards causing injury.

Danger!

For thigh amputees you must use anti-tip tubes.

Danger!

Before setting off, check that your tyre pressure is correct. For rear wheels it should be at least 3.5 bar (350 kPa). The max. pressure is indicated on the tyre. The knee-lever brakes will only work if there is sufficient tyre pressure and if the correct setting has been made (see the Chapter on "Brakes").

Danger!

If the seat and back sling are damaged, you must replace them immediately.

Danger!

Be careful with fire, in particular with burning cigarettes. Seat and back slings can be set alight.

Danger!

If the wheelchair is subject to direct sunlight for a long period of time, then parts of the wheelchair (e.g. frame, legrests, brakes and sideguard) may become hot (>41°C).

Danger!

Always make sure that the quick-release axles on the rear wheels are set properly and lock in. If the button on the quick-release axle is not pressed in, the rear wheel cannot be removed.

Note !

The effect of the knee-lever brake as well as the general driving characteristics are dependant on tyre pressure. The wheelchair is significantly lighter and easier to manoeuvre when the rear wheels are pumped up correctly and both wheels have the same pressure.

Note !

Make sure that your tyres have sufficient tread! Please note that you are subject to any and all traffic laws when driving in public traffic.

Note !

Always be careful with your fingers when working or adjusting the wheelchair!

The products shown and described in this manual may not be exactly the same in every detail as your own model. However, all instructions are still entirely relevant, irrespective of detail differences.

The manufacturer reserves the right to alter without notice any weights, measurements or other technical data shown in this manual. All figures, measurements and capacities shown in this manual are approximate and do not constitute specifications.

Guarantee

Guarantee

THIS DOES NOT AFFECT YOUR LEGAL RIGHTS IN ANY WAY.

Guarantee conditions

1) Repair or replacement is carried out by the authorised Sunrise Medical dealer.

2) To fulfil the guarantee conditions, should servicing need to be carried out on your wheelchair under this agreement, contact the designated Sunrise Medical Authorised Service Agent immediately, with precise details on the type of difficulty. Should you be using the wheelchair outside the area covered by the designated Sunrise Medical customer service agent, the work will be carried out under "guarantee conditions" by another agency as designated by the manufacturer.

3) Should a part or parts of the wheelchair require repair or replacement within 24 months (5 years for frame and cross-brace) after transfer of ownership to the original purchaser, and provided that this person is still the owner of the wheelchair, as a result of a specific manufacturing or material defect, the part or parts will be repaired or replaced free of charge if the wheelchair is returned to the authorised customer service agent.

Note: This guarantee cannot be transferred.

4) The guarantee also covers all repaired or replaced parts for the remaining period of the guarantee for the wheelchair.

5) For spare parts which are fitted after the start of the original guarantee, we give a further 24 month guarantee.

6) Consumable parts are normally excluded from the guarantee, except for any premature wear of these parts which was caused by the original manufacturing fault. These parts include, amongst others, upholstery, tyres, inner tubes and similar parts.

7) The guarantee conditions above cover all product parts for models which were purchased at full sales price.

8) Normally we do not accept responsibility if a repair or replacement of the wheelchair is required for one of the following reasons:

a) The product or part has not been maintained or serviced in accordance with the manufacturer's recommendations as shown in the User Instructions and/or the Service Instructions. Accessories have been used which are not specified as original accessories.

b) The wheelchair or a part of the wheelchair was damaged through neglect, accident or improper use.

c) Alterations to the wheelchair or parts, which are not in accordance with the manufacturer's specifications or the carrying out of repairs before informing the customer service agent.

Wheelchair components

We at SUNRISE MEDICAL have been awarded the ISO 9001 certificate, which affirms the quality of our products at every stage, from R & D to production.

This product complies with the standards set forth in EU directives. Optional equipment and accessories are available at extra charge.

Wheelchairs:

1. Push handles
2. Backrest upholstery
3. Sideguard
4. Seat sling
5. Footrest
6. Castors
7. Footboard
8. Fork
9. Quick-release axle
10. Wheel locks
11. Handrim
12. Rear wheel

Handling

Quick-release axles on rear wheel

The rear wheels are equipped with quick-release axles. The wheels can thus be installed or removed without using tools. To remove a wheel, simply depress the quick-release button on the axle (1) and pull it out (Fig.3).

CAUTION!
Hold the quick-release button on the axle

depressed when inserting the axle into the frame to mount the rear wheels.

Release the button to lock the wheel in place. The quick-release button should snap back to its original position.

Transporting the chair

Removing the rear wheels will keep the chair as compact as possible.

The backrest can be folded down by pulling the cord located on the backrest.

Options – Step Tubes

Step Tubes

Step tubes are used by attendants to push a wheelchair over an obstacle. Simply step on the tube to push a wheelchair, for example, over a kerb or step (Fig. 4).

Options - Brakes


CAUTION – Braking power can be affected by incorrect fitting and adjustment of the brakes, as well as tyre pressure which is too low.

Wheel locks

Your wheelchair is equipped with two wheel locks. They are applied directly against the tyres. To engage, press both wheel lock levers forward against the stops. To release the wheel lock, pull the levers back to their original positions. Braking power will decrease with:

- Worn tyre tread
- Tyre pressure that is too low
- Wet tyres
- Improperly adjusted wheel locks.

The wheel locks have not been designed to be used as brakes for a moving wheelchair. The wheel locks should therefore never be used to brake a moving wheelchair. Always use the handrims for braking. Make sure that the interval between the tyres and wheel locks complies with the given specifications. To readjust, loosen screw (1) and set the appropriate interval. Tighten screws again (see the page on torque) (Fig.5).

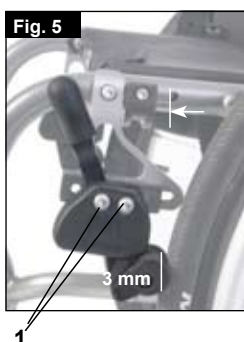Standard Wheel lock (Fig.6)

Compact brakes

Compact brakes are underneath the seat sling and are operated by pulling the brakes towards the rear, in the direction of the tyre. For the brakes to work properly, this must be pulled until it reaches the stops (Fig.7).

Safari wheel lock (Wheel lock integrated into Side guard):

The Safari wheel lock are operated by pushing the lever forwards (away from the patient's body). To release the brakes, pull the lever back towards the body (Fig.8).

Brake lever extension

The longer lever helps to minimise the effort needed to set the wheel locks. The brake lever extension is screwed to the brakes. By raising this, it can be flipped forward (Fig.9).

CAUTION!
Incorrect mounting of the wheel lock will result in a higher effort to operate. This might cause the wheel lock extension lever to break.

Options – Footplate Adjustment

CAUTION – Do not stand on the footboard! Even if the user is sitting in the chair, there is still a risk of tipping over and injury.
CAUTION – When transferring, do not stand on the footboard, there is a risk of tipping over and injury.

Adjusting the footrest

Removing the screw (1) will allow you to adjust the footrest to correspond to the length of your lower leg and re-attach the footrest. The angle of the footrest may be individually adjusted (not for footrest with bar) by loosening screws (2).

The side protection on the footrest prevents the feet slipping off accidentally. Make sure that after any adjustment work, all screws are tightened correctly (see the page on torque) (Fig.10).

Flip-up footrest:

The footrest can be flipped up to the rear and this can make it easier for the patient to transfer to/from the chair. The footrest is flipped up by raising the front section. To return the footrest to its original position, the front section must be flipped downwards (Fig.11).

High-mount footrest:

The high-mount footrest is fitted on the inner part of the frame and permits a higher footrest position (Fig.12).

Options - Castor

Castor, Castor adapter, Castor fork

From time to time the wheelchair may veer slightly to the right or left or the castors may flutter. This may be caused by the following:

- The forward or reverse wheel motion has not been set properly.
- The castor angle has not been adjusted properly.
- Castor and/or rear wheel air pressure is incorrect; wheels do not turn smoothly.


The optimum adjustment of the castors is required for the wheelchair to run in a straight line.

Castors should always be adjusted by an authorised dealer. The castor plates must be readjusted, and the wheel locks must be checked any time the rear wheel position has been altered.

Adjusting the castor


To ensure that both forks are set parallel, simply count the teeth visible on both sides. After setting the castor fork, the teeth will guarantee a secure position, allowing an adjustment of 16° in 2° increments. Use the flat side to check for a right-angled position to the ground.

The patented design allows the castor fork to be turned, so that it can be reset at right-angles to the ground when the seat angle is changed (Fig. 13 + 14).

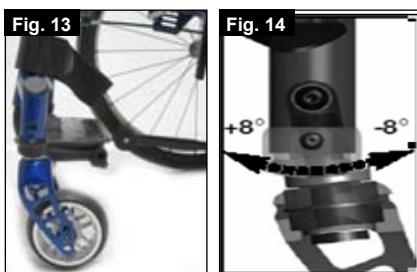Setting the directional stability

Remove the cover cap (1). Then release the Allen screws (2) on the inside. You can then undo the screw (3). You can now turn the black socket (4) left or right. Left – chair pulls to the left Right – chair pulls to the right Then tighten up the screws (3) again. Please set a 90° angle from the fork to the floor. Now tighten up the screws (2) again and close the tube with the cover cap (1) (Fig.15)

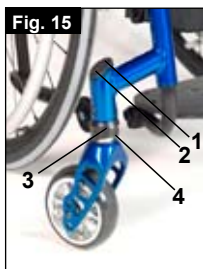Options - Seat

Adjusting the seat height

To adjust the rear seat height release and remove the Allen screw (1) (one on each side), which fix the clamp of the camber tubes (2) to the axle plate (3). Adjust the two camber tube clamps (2) to the required height and replace the two Allen screws. Before tightening the screws, please follow the instructions for setting the toe-in/toe-out to zero.

Tighten the screws to 10 Nm (Fig.16)

NOTE – An adjustment to the castor angle may be necessary when adjusting the rear seat height.

Options – Wheel Alignment

Adjusting the wheel alignment

Important: To achieve the very best movement, the rear wheels must be adjusted to their optimum position, which means correctly adjusting the wheel alignment.

To do this, measure the distance between both wheels front and rear to ensure that they are parallel to one another. The difference between both measurements should not exceed 5 mm.

To adjust the wheels to make them parallel, loosen the screws and turn the axle sleeve accordingly. Make sure that after any adjustment work, all screws are tightened correctly (see the page on torque).

SIMBA tracking adjustment


Setting the toe-in/toe-out to zero

NOTE: A wheelchair with 0° camber cylinders cannot have toe-in or toe-out. This setting is necessary only with 3°, 6° and 9° camber cylinders.

The term "toe-in or toe-out" defines how well the rear wheels of the chair are aligned in relation to the ground. This determines how well the chair will run. Normal resistance or rolling resistance is present when toe-in is set to zero.

To set toe-in/toe-out to zero: Loosen the Allen screws (1) (one on each side), which secure the angle tube clamp. Check the ball in the horizontal (2) plane and turn the angle tube (3) until the ball is in the centre. Toe-in is now zero.

Before tightening the screws (1), check that the camber tube is centred left-to-right. The gap should be the same on both sides, or there should be no gap at all. Tighten the screws to 10 Nm (Fig.17 + 18 + 19).

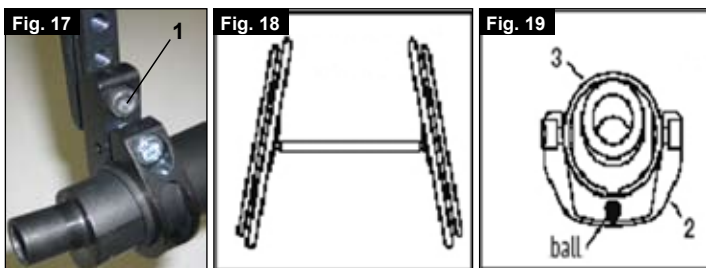Adjusting the wheelbase width


K. REAR WHEELBASE

The rear wheelbase is defined as the distance between the upper side of the rear wheels and the backrest tubes, and is represented by measurement X. The factory setting is (1.25 cm). This has to be increased if a larger gap between the tyres and the optional height-adjustable armrests has to be created.

NOTE: When adjusting the rear wheelbase, adjust first one wheel then the other. If both sides are loosened at the same time, this will alter the toe-in/toe-out adjustment. To adjust the rear wheelbase, the parts of the camber (4) move telescopically into or out of the camber tube (5), and lock into place when they reach the end. Loosen screw (6) (located closest to the camber tube) on the left side of the chair.


Move the quick-release axle inwards or outwards to achieve the desired wheelbase. Tighten the screws to 10 Nm. Repeat this procedure on the right side of the chair and adjust the gap so that it is the same amount as on the left side (Fig.20).

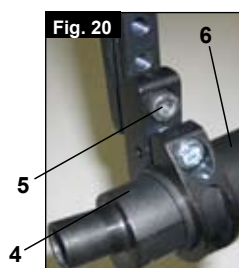Options – Backrest

Angle-adjustable back

The backrest angle is adjusted by altering the position of the backrest mounting. The screws (1) must be loosened and removed on both sides. The backrest mounting should now be adjusted until the desired angle is reached. Then the screws (1) must be re-tightened and secured (Fig.21)

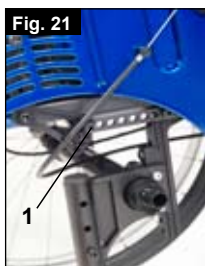CAUTION – The screws must be re-tightened. Otherwise the angle adjustment will be lost.

Folding backrest

Release the backrest by pulling the cord and, at the same time, push it forward to fold it down. To return the backrest to its original position, this must be pulled back as far as possible until it locks into place on both sides.

NOTE – If the backrest is pulled too far backwards, operating the brakes could cause the user to fall backwards. Therefore, please make sure that both backrests lock into place.
CAUTION – When folding the backrest down, please make sure that your fingers do not get caught.

Backrest in the relax position

A relax position can be achieved with a 10° backrest angle. Pull the cord located on the backrest to release it, and pull it backwards until both backrest tubes lock into place.

NOTE – If the backrest is pulled too far backwards, operating the brakes could cause the user to fall backwards. Therefore, please make sure that both backrests lock into place.
CAUTION – When folding the backrest down, please make sure that your fingers do not get caught.

Adjustable back sling

The adjustable back sling can be adjusted for tension by using several straps. The back sling upholstery can be accessed from the inside via an opening and can be padded to suit individual tastes.

Height-adjustable backrest

The back height may be set to various back heights, in 2.5 cm steps. The adjustment ranges are 20-25 cm, 25-30 cm, 30-35 cm and 35-40 cm. Release the screws (1) and set the backrest to the desired height. Tighten up the screws again (see the page on torque) (Fig.22).

Options - Sideguard

The sideguards of the SIMBA are primarily designed for the new integrated sideguard brakes. If you have not ordered the integrated sideguard brakes then the openings in the sideguard are fitted with covers.

CAUTION – If the sideguard is used without the covers, then this presents a risk of injury. Where the integrated sideguard brake is not fitted, then there must be a cover.

Armrest

The armrest is screwed into place on the sideguard at the designated position. The armrest is height-adjustable.

Therapy table

A therapy table can only be used in combination with the armrest. It is fixed onto the armrest and is positioned on both armrests.

Options – Push handles

CAUTION – A wheelchair should never be lifted with an occupant; it should only be pushed.

Height-adjustable push handles

These handles are secured with pins to prevent them from sliding out unintentionally. Opening the quick-release lever (1) makes it possible to adjust the height of the push handles (2) to meet your individual needs. As you move the lever, you will hear a locking mechanism; you may now easily position the push handle as desired.

The nut on the tension lever determines how tightly the push handles are clamped into place. If the nut is loose after adjusting the tension lever, the push handle will also be too loose. Turn the push handle from side to side before use to make sure that it is clamped securely enough into place.

After adjusting handle height, always clamp the tension lever (1) securely into place. If the lever is not secure, injuries could result when ascending stairs (Fig. 23 + 24).


NOTE – If the height-adjustable push handles are not fitted properly, there is a risk that these will develop "play" or that they move out of position. Please make sure that the relevant screws are tightened correctly.

Fold-down push handles

If the push handles are not in use, they can be folded down depressing the button (2). When they are needed again, simply flip them back up until they click into place (Fig.25).

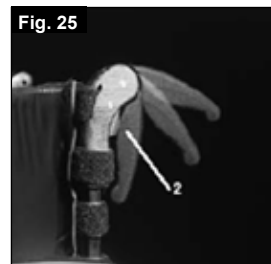Options – Positioning belt


The positioning belt is a hip belt that enhances wheelchair safety. The positioning belt is fastened to the frame as shown in the picture.

CAUTION - This positioning belt must not be used as a safety belt in a car!

Options – Anti-Tip Tubes

Sunrise Medical recommends anti-tip tubes for all chairs. When fitting anti-tip tubes, use a torque of 16 Nm (Fig.26)

1. Slotting the anti-tip tubes into the clamp.
 - a. Press the rear button on the anti-tip tube on the anti-tip tube adaptor, so that both release pins are drawn inwards.
 - b. Slot the anti-tip tubes (1) into the anti-tip tube adaptor.
 - c. Turn the anti-tip tubes downwards until the release pin locks into the clamp.
 - d. Fit the second anti-tip tube in the same way.

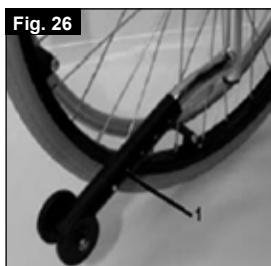2. Adjusting the anti-tip tubes
To achieve the correct ground clearance of approx. 1" to 2" (2.5 cm to 5.0 cm), the anti-tip tubes must be raised or lowered. Press the anti-tip tube release button, so that both release pins are drawn inwards. Move the inner tube up or down to slot into the height holes provided. Release the button. Fit the second anti-tip wheel in the same way. Both wheels should be at the same height.

Danger!

If the anti-tip tubes are not fitted, or have been fitted incorrectly, there is a risk of tipping over and of injury.

Options - Seat

Seat sling

Remove the screws on the left side of the sling. Adjust the VELCRO® material to increase the tension of the seat sling. Tighten up the screws again. If the screws are difficult to fit back in place, try and locate the holes with a sharp object. Ensure also that the plastic base is in the correct position before the screws are tightened again.

Options – Seat width and seat depth growth

Seat depth growth

The seat width can be expanded by one step. To achieve this:

- The screws in the footrest and axle, as well as the two stabilising bars, must be loosened,
- Back and seat slings must be loosened using the Velcro and set to the next largest seat width,
- Footrest, axle and the two stabilising bars must be set to the next largest seat width and then re-tightened.

Seat depth growth

The seat depth can be increased by two steps. To achieve this:

- Two screws must be loosened and the complete part must be pushed backwards by one or two positions,
- The screws must then be re-tightened.

Options – Crutch holder

Crutch holder

This device permits crutches to be transported directly on the wheelchair. It has a Velcro loop to fasten crutches or other aids.

CAUTION!

Never try to use or even remove the crutches or other aids while moving.

Tyres and Mounting

Solid rubber tyres are standard.

With pneumatic tyres always make sure that the tyres have the correct air pressure, as otherwise the performance of the wheelchair may be affected. If the tyre pressure is too low, rolling resistance will increase, requiring more effort to move the chair forward. Low tyre pressure also has a negative impact on manoeuvrability. If the tyre pressure is too high, the tyre could burst. The correct pressure for a given tyre is printed on the surface of the tyre itself.

Tyres are mounted in the same way as an ordinary bicycle tyre. Before installing a new inner tube, you should always make sure that the base of the rim and the interior of the tyre are free of foreign objects. Check the pressure after mounting or repairing a tyre. It is critical to your safety and to the wheelchair's performance that regulation air pressure be maintained and that tyres be in good condition.

Maintenance and care

- Check the tyre pressure every four weeks. Check all tyres for wear and damage.
- Check the brakes approximately every four weeks to make sure that they are working properly and are easy to use.
- Change tyres as you would an ordinary bicycle tyre.
- All of the joints that are critical to using your wheelchair safely are self-locking nuts. Please check every three months to make sure that all bolts are secure (see the section on torque). Safety nuts should only be used once and should be replaced after use.
- Please use only mild household cleansers when your wheelchair is dirty. Use only soap and water when cleaning the seat upholstery.
- If your wheelchair should ever get wet, please dry it afterwards.
- A small amount of sewing-machine oil should be applied to quick-release axles approximately every eight weeks. Depending on the frequency and type of use, we recommend taking your wheelchair to your authorised dealer every six months to have it inspected by trained personnel.

CAUTION:

Sand and sea water (or salt in the winter) can damage the bearings of the front and rear wheels. Clean the wheelchair thoroughly after exposure.

The following parts can be removed and sent back to the manufacturer/dealer for repair:

- Rear wheels
- Armrest
- Footrest holder
- Anti-Tip tubes

These components are available as spare parts. For further information, please see the spare parts catalogue.

Hygiene measures when being re-used:

Prior to the wheelchair being re-used, it must be carefully prepared. All surfaces which come into contact with the user must be treated with a disinfection spray.

To do this, you must use a rapid alcohol-based disinfectant for medical products and medical devices. We recommend AZO Max active wipes.

Please take into account the manufacturer's instructions for the disinfectant you are using.

In general, a complete disinfection cannot be guaranteed on seams. We therefore recommend that you dispose of seat and back slings to avoid microbacterial contamination with active agents according to § 6 infection protection law.

Options – Lap belt

Before using your wheelchair ensure the lap belt is worn.

The lap belt must be checked on a daily basis to ensure it is free from any obstruction or adverse wear. Always make sure that the lap belt is correctly secured and adjusted prior to use. Too loose a strap could cause the user to slip down and risk suffocation or cause serious injury.

The lap belt is fitted to the wheelchair as shown in the illustrations. The seat belt comprises 2 halves. They are fitted using the existing seat stay retaining bolt fitted through the eyelet on the belt. The belt is routed under the rear of the side panel. (Fig. 27)

Adjust the belt position so buckles are in the centre of the seat. (Fig. 28)


Adjust lap belt to suit the user's needs as follows:

To reduce the belt length	To increase the belt length
	
Feed free belt back through male buckle and slide adjusters.	Feed free belt through slide adjusters and male buckle to provide more belt length.

When fastened check space between the lap belt and user, when correctly adjusted it should not be possible to insert more than the flat of the hand between the lap belt and the user. (Fig. 29)

The lap belt should be fixed so that the belt sits at an angle of 45 degrees across the users pelvis. The user should be upright and be as far back as possible in the seat when correctly adjusted. The lap belt should not allow the user to slip down in the seat. (Fig. 30)

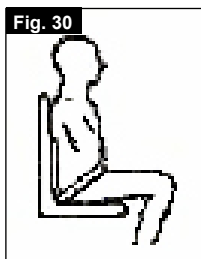To fasten buckle:

Firmly push male buckle into female buckle.

To release belt:

Press exposed sides of male buckle and push towards centre whilst gently pulling apart.

If in doubt about the use and operation of the lap belt then ask your healthcare professional, wheelchair dealer, carer or attendant for assistance.

Do not rely on the lap belt only when the wheelchair and occupant are transported in a vehicle. Use the separate occupant lap and diagonal restraints provided in the vehicle. See section transportation for more details.

Advice to client

The lap belt must only be fitted by an approved Sunrise Medical dealer / agent. The lap belt should only be adjusted by a professional, or a Sunrise Medical approved dealer / agent. The lap belt must be checked on a daily basis to ensure they are adjusted correctly and are free from any obstruction or adverse wear.

Sunrise Medical does not encourage the transportation of any person in a vehicle using this lap belt as a method of restraint. Please see Sunrise Medical transit booklet for further advice on transportation.

Maintenance:


Check lap belt, and securing components, at regular intervals for any sign of frays, or damage. Replace if necessary.

NOTE:

The lap belt should be adjusted to suit the end user as detailed above. Sunrise Medical recommend that the length and fit of the belt be checked on a regular basis to reduce the risk of the end user inadvertently re-adjusting the belt to an excessive length.

Disposal / Recycling of materials

If the wheelchair has been made available to you free of charge, then it does not belong to you. If it is no longer required, then follow the instructions to return it as given by the organisation that made the wheelchair available to you.

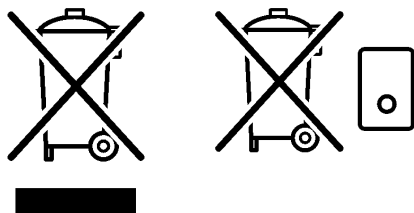In the following section, there is a description of the materials used on the wheelchair, in view of the disposal or recycling of the wheelchair and its packaging.

Particular regulations with regard to disposal or recycling may be in force locally and these must be taken into account when performing disposal. (This can include the cleaning or decontamination of the wheelchair prior to disposal).

Aluminium: Castor forks, wheels, sideguards for the chassis, armrest frame, footrest, push handles

Steel: Fixing points, quick-release axle

Plastic: Handles, tube stoppers, castors, footplates, armpads and 12" wheel/tyre

Packaging: Plastic bags made of soft polyethylene, cardboard

Upholstery: Polyester fabric with PVC coating and highly-flammable foam.

Disposal or recycling should be carried out by a disposal company or at a public disposal point. You can also return your wheelchair to your dealer for disposal.

Trouble-shooting

Wheelchair pulls to one side

- Check tyre pressure
- Check to make sure wheel turns easily (bearings, axle)
- Check the castor angle
- Check to make sure both castors are making proper contact with the ground

Castors begin to wobble

- Check the castor angle
- Check to make sure all bolts are secure; tighten if necessary (see the section on torque)
- Check to make sure both castors are making proper contact with the ground

Wheelchair / Cross-tube assembly does not snap into position in the seat saddle

- Chair is still new, i.e., the seat or backrest upholstery is still very stiff. This will improve with time.

Wheelchair is difficult to fold up

- Adjustable backrest upholstery is too stiff. Loosen it accordingly.

Wheelchair squeaks and rattles

- Check to make sure all bolts are secure; tighten if necessary (see the section on torque)
- Apply small amount of lubrication to spots where movable parts come in contact with one another

Wheelchair begins to wobble

- Check angle at which castors are set
- Check tyre pressure
- Check to see if rear wheels are adjusted differently

Nameplate

Nameplate

The nameplate is located on either the cross-tube assembly or the transverse frame tube, as well as on a label in the owner's manual. The nameplate indicates the exact model designation and other technical specifications. Please provide the following pieces of information whenever you have to order replacement parts or to file a claim:

- Serial number
- Order number
- Month/Year

Technical Data

Overall width:

With 20" standard wheels, incl. handrims with 6° camber: SIMBA: SW + 29 cm

Overall length: 88.9 cm with seat depth 36

Overall height: 87.7 cm with backrest height 40

Weight in kg: from 9 kg






Maximum load:

SIMBA up to a load of 65 kg

Seat heights:

The choice of frames, forks and castors, as well as the rear wheel size (20", 22", 24") determines the seat heights which can be achieved.

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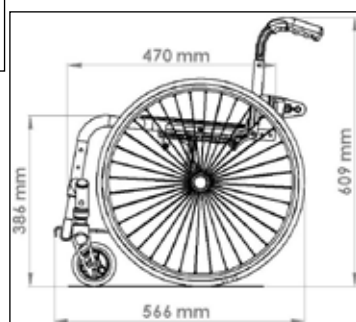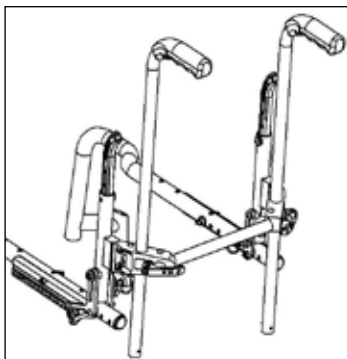
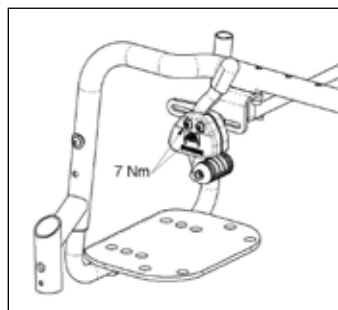
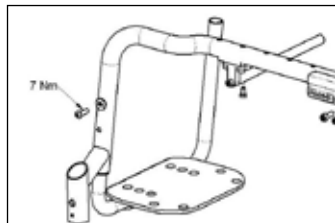
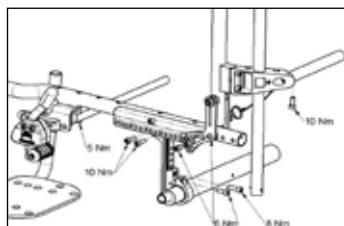
Technical drawing of a mechanical assembly. The drawing shows a horizontal beam with a vertical support on the left and a horizontal support on the right. A vertical rod is attached to the left end of the beam, and a horizontal rod is attached to the right end. Dimensions are indicated: 10 Nm for the vertical rod, 10 Nm for the horizontal rod, 10 Nm for the vertical support, 10 Nm for the horizontal support, 10 Nm for the vertical rod, and 10 Nm for the horizontal rod.

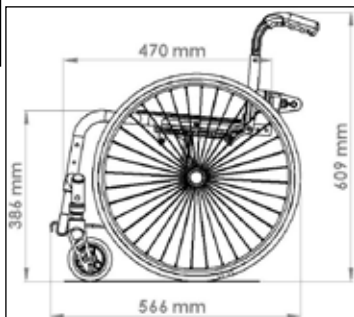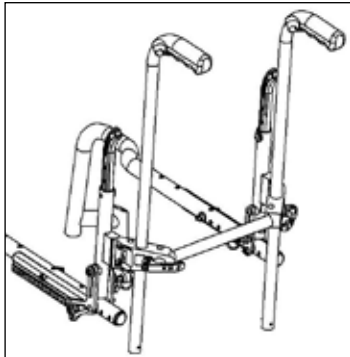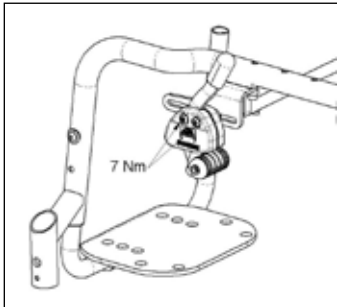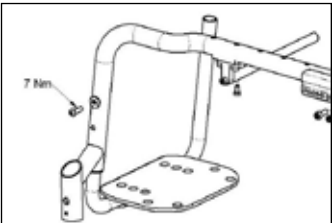Baghjul	Bageste sæde-højde
20"	40
20"	38
20"	37
20"	36
20"	35
20"	34
20"	32
20"	31
20"	30
20"	29
22"	43
22"	42
22"	41
22"	40
22"	38
22"	37
22"	36
22"	35
22"	34
22"	33
24"	46
24"	45
24"	43
24"	42
24"	41
24"	40
24"	39
24"	38
24"	36
24"	35

Tekniske data

DANSK

KORT stelversion (230 mm)			LANG stelversion (270 mm)		
Svinghjul	Forgaffel 45 mm bred	Forreste sædehøjde i cm	Svinghjul	Forgaffel 45 mm bred	Forreste sædehøjde i cm
3"	111 mm	37	3"	111 mm	41
3"	111 mm	38	3"	111 mm	42
3"	111 mm	40	3"	111 mm	44
3"	123 mm	38	3"	123 mm	42
3"	123 mm	39	3"	123 mm	43
3"	123 mm	40	3"	123 mm	44
4"	111 mm	37	4"	111 mm	41
4"	111 mm	38	4"	111 mm	42
4"	111 mm	39	4"	111 mm	43
4"	111 mm	40	4"	111 mm	44
4"	123 mm	38	4"	123 mm	42
4"	123 mm	39	4"	123 mm	43
4"	123 mm	40	4"	123 mm	44
4"	123 mm	41	4"	123 mm	45
4"	123 mm	42	4"	123 mm	46
5"	111 mm	39	5"	111 mm	43
5"	111 mm	40	5"	111 mm	44
5"	111 mm	41	5"	111 mm	45
5"	123 mm	41	5"	123 mm	45
5"	123 mm	42	5"	123 mm	46
5"	123 mm	43	5"	123 mm	47
6"	123 mm	43	6"	123 mm	47
6"	123 mm	44	6"	123 mm	48
KORT stelversion (230 mm)			LANG stelversion (270 mm)		
Svinghjul	Forgaffel 32 mm bred	Forreste sædehøjde i cm	Svinghjul	Forgaffel 32 mm bred	Forreste sædehøjde i cm
3"	98 mm	36	3"	98 mm	40
3"	98 mm	37	3"	98 mm	41
3"	98 mm	38	3"	98 mm	42
3"	111 mm	37	3"	111 mm	41
3"	111 mm	38	3"	111 mm	42
3"	111 mm	39	3"	111 mm	43
4"	98 mm	37	4"	98 mm	41
4"	98 mm	38	4"	98 mm	42
4"	98 mm	39	4"	98 mm	43
4"	111 mm	38	4"	111 mm	42
4"	111 mm	39	4"	111 mm	43
4"	111 mm	40	4"	111 mm	44
5"	98 mm	40	5"	98 mm	44
5"	98 mm	41	5"	98 mm	45
5"	111 mm	40	5"	111 mm	44
5"	111 mm	41	5"	111 mm	45
5"	111 mm	42	5"	111 mm	46

Sunrise Medical GmbH & Co.KG
Kahlbachring 2-4
69254 Malsch/Heidelberg
Germany
Tel.: +49 (0) 7253/980-0
Fax: +49 (0) 7253/980-111
www.sunrisemedical.com

Sunrise Medical Limited
High Street Wollaston
West Midlands DY8 4PS
England
Tel.: +44 (0) 1384446688
www.sunrisemedical.com

Sunrise Medical S.L.
Polígono Bakiola, 41
48498 Arrankudiaga – Vizcaya
España
Tel.: +34 (0) 902142434
Fax: +34 (0) 946481575
www.sunrisemedical.com

Sunrise Medical S.A.
13, Rue de la Painguetterie
37390 Chauceaux sur Choissille
France
Tel.: +33 (0) 247554400
Fax: +33 (0) 247554403
www.sunrisemedical.com

Sunrise Medical Srl
Via Riva 20, Montale
29100 Piacenza (PC)
Italy
Tel.: +39 0523-573111
Fax: +39 0523-570060
www.sunrisemedical.com

Sunrise Medical A.G.
Lückhalde 14
3074 Muri bei Bern
Switzerland
Tel.: +41 (0) 31-958-3838
Fax: +41 (0) 31-958-3848
www.sunrisemedical.com

Sunrise Medical AS
Rehabsenteret
1450 Nesoddtangen
Norway
Tel.: +47 (0) 66963800
Faks: +47 (0) 66963838
www.sunrisemedical.com

Sunrise Medical AB
Britta Sahlgrens gata 8A
421 31 Västra Frölunda
Sweden
Tel: +46 (0) 31-748 37 00
Fax: +46 (0) 31-748 37 37
www.sunrisemedical.com

Sunrise Medical B.V.
Groningenhaven 18-20
3433 PE Nieuwegein
Nederlands
Tel.: +31 (0) 3060-82100
Fax: +31 (0) 3060-55880
www.sunrisemedical.com

G2 A/S
Graham Bells Vej 21-23 A
8200 Aarhus N
Denmark
Tel.: +45 87 41 31 00
Fax: +45 87 41 31 31
www.G2.dk